## Körperanalysewaage Tabea/Feline

Body Analyser Scale | Báscula de análisis corporal

Balance d'analyse corporelle | Bilancia con analisi corporea

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Deutsch

English

Español

Français

Italiano

Polski

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BA 1300 | BA 1301 | BA 1302 | BA 1303 | BA 1304

Operating instructions | Instrucciones de operación | Mode d'emploi

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# Dear Customer!

You have chosen to purchase a high-quality product from the brand ADE, which combines intelligent functions with an exceptional design.

This set of body analyser scales enables you to measure and get to know yourself and your body better.

The many years of experience of the brand ADE ensures a high technical standard and proven quality.

We wish you much enjoyment while weighing yourself and discovering your values!

Your ADE Team

## General Information

---

### About this manual

This operating manual describes how to safely use and care for your scale. Keep this operating manual in a safe place in case you would like to refer to it later. If you pass these scale on to someone else, pass on this manual as well.

Failure to observe the recommendations in this operating manual can lead to injuries or damage to your scale.

### Explanation of symbols

This symbol, combined with the word DANGER, warns against serious injury.

This symbol, combined with the word WARNING, warns against moderate and mild injuries.

This symbol, combined with the word NOTE, warns against material damage.

This symbol signifies additional information and general advice.

# Table of Contents

---

General Information.....	16
Table of Contents .....	17
Intended Use.....	18
Safety .....	18
Overview .....	20
Scope of delivery .....	20
Start-up .....	21
This is how your body analyser scale work .....	21
Measuring your weight .....	22
Entering personal data.....	23
Body analysis.....	23
Changing the battery .....	26
Cleaning.....	26
Fault/Remedy.....	27
Technical Data .....	28
Declaration of Conformity .....	28
Guarantee .....	29
Disposal .....	29

## Intended Use

---

The bathroom scale is designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

Likewise, the body analysis function is designed only for private use. The scale is not suitable for use in the professional sporting sphere.

The maximum carrying capacity of 180 kg must not be exceeded.

Do not use the scale if you ...

- ... are pregnant,
- ... have a cardiac pacemaker,
- ... are a dialysis patient,
- ... suffer from oedema,
- ... have implants, e.g. artificial limbs, metal plates, etc.

**If in doubt, consult your GP.**

## Safety

---

- This item can be used by children 8 years of age and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge, if they are supervised or have been instructed in the safe use of the item and understand the dangers that may arise.
- Children must not play with the item.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries can be life-threatening. Therefore, keep batteries out of reach of small children. If a battery has been swallowed or has entered the body in any other way, seek medical help immediately. Failure to do so within 2 hours may result in severe internal burns that could cause death.

## DANGER for children

- Keep children away from packaging material. Risk of choking if swallowed. Also make sure that children do not pull the packaging material over their heads.

## **RISK of fire/burning and/or explosion**

- There is a risk of explosion if batteries are replaced incorrectly. Replace the battery using only the same or an equivalent type of battery. Pay attention to the "Technical Data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited.
- Avoid contact with skin, eyes, and mucosa. If you come into contact with battery acid, flush the affected site immediately with plenty of clean water and see a doctor promptly.

## **WARNING against risk of injury**


- Do not step onto the outer edge of any one side of the scale. Risk of tipping!
- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!
- The weighing surface is made of glass. If you fall down, there is a risk of the glass breaking and thus a risk of injury from glass splinters.

## **BEWARE of material damage**

- Place the scale on a stable, even surface.
- Do not overload the scale. Pay attention to the "Technical Data". Place yourself on the scale exclusively for weighing purposes.
- Protect the scale from sharp jolts and vibrations.
- Remove the battery from the scale once it has been used or when you are no longer using the scale. You will thus avoid damage that may result from the leakage of battery acid.
- Pay attention to the correct polarity (+/-) when exchanging the battery.
- Do not subject the battery to any extreme conditions, by – for example – storing it on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the scale. Only allow repairs to be carried out at a specialist workshop. Repairs carried out incorrectly can considerably endanger the user.
- There are feet on the underside of the scale. Since floor surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place an anti-slip underlay beneath the scale.

# Overview

---

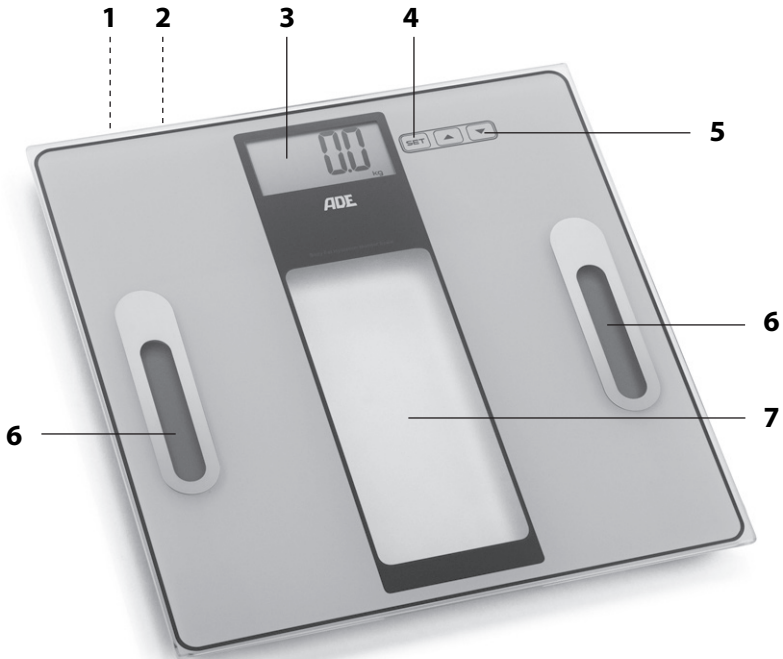- 1 Battery compartment (reverse side)
- 2 **kg lb st** button (reverse side): Changing the unit of measurement
- 3 Display
- 4 **SET** button: Call up setting mode, store entries, call up stored measurement data
- 5 Buttons **▲ ▼**: Select memory cell, reduce data, increase data
- 6 Sensor
- 7 Weighing surface

## Scope of delivery

---

- Scale
- Battery Type CR2032, 3 V  (already inserted)
- Operating instructions

## Symbols on the scale

Auto-on

The scale automatically switch on when you step onto the weighing surface with a little force.

CAUTION:

Slippery when wet !

**WARNING:** When the weighing surface is damp, there is a risk of slipping.

## Start-up

### Removing the iso-strips

When the scale is delivered, the battery is already inserted and fitted with iso-strips to prevent it from discharging ahead of time. In order to start up the scale, you must remove the iso-strips.

The battery compartment is located on the underside of the device.

1. With a little force, press on the arrow on the battery compartment cover and push it upwards and off.
2. Pull out the iso-strips.
3. Push the battery compartment cover back into place.

### Choosing the unit of measurement

When the scale is delivered, the unit of measurement is set to “kg” as standard. However, the measurement result can also be displayed in the English units of measurement, “lb” (pound) or “st:lb” (stone:pound).

- Press the **kg lb st** button next to the battery compartment on the reverse side of the scale repeatedly. Pay attention to the display indicator.

### Positioning the scale

- Place the scale on an even, firm, and dry surface. Do not place the scale on carpet.

## This is how your body analyser scale work

Your new body analyser scale works using so-called BIA-technology (Bioelectric Impedance Analysis) This involves an undetectable, completely harmless and weak current flowing through your body.

The measurement takes place via sensors that are mounted on the weighing surface.

- The current flows relatively unhindered through the fluid contained in muscle mass.
- The resistance of fat mass is stronger. This resistance is measured.

Using your personal data (sex, weight, age, and height), the scale calculates the following information about your body composition:

- Weight
- Body fat proportion in %
- Body water in %
- Muscle mass in %
- Bone mass (kg), and your daily
- Basal metabolism in Kcal (energy demand)

The scale has a total of 12 user slots.

## Measurement errors

Measurement or calculation errors may arise for the following groups of people during the body analysis:

People who ....

... spend 10 hours or more per week playing sport,

... are professional athletes, body builders or endurance athletes,

... are over 70 years old,

... are less than 10 years old.

The “normal” weighing function (see next chapter) is not affected by these considerations.

## Measuring your weight

---

For normal weighing – without body analysis – you do not have to enter any data.

- Position yourself on the scale. Remain standing still and do not hold on to anything.

The measured weight will be displayed shortly afterwards.

A few seconds after you step off the weighing surface, the scale will switch off automatically.

## Entering personal data

1. Briefly step onto the weighing surface once, using a little force.  
"0.0" and the set unit of measurement, e.g. "kg", will be shown on the display.
2. Press SET.  
On the display, a flashing memory user slot (1-12) will be shown.
3. Using the ▲ or the ▼ button, select the desired user slot.
4. Press SET.  
Using the ▲ or the ▼ button, choose your sex: ♀ for male, or ♂ for female.
5. Press SET.  
On the display, the preset value for height (e.g. "165" for cm) will flash.
6. Using the ▲ or the ▼ button, enter your height.
7. Press SET.  
On the display, the preset value for age (preset: 25) will flash.
8. Using the ▲ or the ▼ button, enter your age.
9. Finally, press SET. You have thus entered the required data.  
On the display, a user slot number will flash again and you can proceed, as previously described, and enter the data of another person.  
Or wait a few seconds until "0.0" is shown on the display and you can undertake your body analysis.  
If you do nothing else, the scale will automatically switch off after a few seconds.

- If you hold down the ▲ or the ▼ button, the forward and backward flow of the displayed data will occur more quickly.
- Entered data can only be overwritten by making a new entry.

## Body analysis

- Body analysis is only possible if your data has been entered (see previous chapter).
  - See the chapter "Weighing Tips".
1. Briefly step onto the weighing surface once, using a little force.  
"0.0" and the set unit of measurement, e.g. "kg", will be shown on the display.
  2. Briefly press SET.
  3. Using the ▲ or the ▼ button, select the user slot in which you have recorded your personal body data.  
The data is briefly displayed, then "0.0" and the selected unit of measurement, e.g. "kg".

4. Position yourself barefoot on the scale. Remain standing still and do not hold on to anything.
  - In doing so, ensure each foot is placed on a sensor.
  - There must be no contact between your feet, your legs, your calves, or your thighs. Otherwise, the measurement will be imprecise.

The measured weight will be displayed shortly afterwards. The indicator “o” then runs across the display; the body analysis is being carried out. Remain standing still.

Your measured values will then be shown to you repeatedly on the display:

Symbol	Meaning
1.	Your body weight in the set unit of measurement
2. 	Your body fat proportion in % with simultaneous evaluation indicator:  = too little body fat  = slender  = normal  = too much body fat  = obese
3. 	Your body water in %
4. 	Your muscle mass in %
5. 	Your bone mass in kg
6. 	Your basal metabolism in kcal

## Calling up stored values

After an analysis has been performed, you can display the values of your last measurements at any time.

1. Briefly step onto the weighing surface once, using a little force.  
“0.0” and the set unit of measurement, e.g. “kg”, will be shown on the display.
2. Press SET and select your memory cell.  
The memory cell indicator stops flashing after a few seconds and the indicator “0.0” appears on the display, as well as the unit of measurement, e.g. “kg”.
3. By pressing ▲ or ▼ once, select the automatic display or, by repeatedly pressing ▲ or ▼, select the manual display of the last set of measurements.
4. Press SET in order to switch the scale back into analysis mode.

## Weighing tips

- Weighing for the body analysis can only be carried out barefoot. In doing so, your feet should only be lightly moistened. Completely dry feet can lead to false results, as the electrical conductivity will be poor then.

- During the measuring process, stand still and do not hold on to anything.
- Do not weigh yourself directly after unusual physical exertion.
- Do not weigh yourself immediately after standing up. Wait approximately fifteen minutes. In this time, the water is able to distribute through your body.
- Always attempt to weigh yourself at the same time, in order to obtain comparable values.

## Notes on the ascertained values

The proportion of body fat gives the amount of fat as a percentage of the total weight of the body. A universal recommendation for this value does not exist, since it depends significantly on sex (because of different physiques) and age.

Body fat is vital for day-to-day body function. It protects organs, cushions joints, regulates body temperature, stores vitamins, and serves the body as an energy store.

For women, a body fat proportion between 10% and 13% is essential; for men, between 2% and 5% with the exception of competitive sportsmen.

The proportion of body fat is not an unambiguous indicator of health. Excessive body weight and a too-high proportion of body fat are linked with the appearance of many diseases of civilisation such as diabetes, cardiovascular illnesses, etc.

Physical inactivity and poor diet are often triggers for these diseases. There is a clear causal connection between excess weight and physical inactivity.

A considerable part of the human body consists of water. This body water is variably distributed throughout the whole body. Lean muscle tissue contains about 75% water, blood around 83% water, body fat approx. 25% water and bones have a water proportion of about 22%.

- In men approx. 60% of body mass consists of water.
- In women, that value is about 55% (determined by a higher proportion of body fat).

Of this total body water approx. two thirds is stored in the cells and therefore called intracellular water. The other third is extracellular water.

Women 					
Age	Body fat				Body water
	Slender	Normal	Slightly overweight	Overweight	Normal
10-16	< 18 %	18-28 %	29-35 %	> 35 %	57-67 %
17-39	< 20 %	20-32 %	33-38 %	> 38 %	47-57 %
40-55	< 23 %	23-35 %	36-41 %	>41 %	42-52 %
56-85	< 24 %	24-36 %	37-42 %	>42 %	37-47 %

Men 					
Age	Body fat				Body water
	Slender	Normal	Slightly overweight	Overweight	Normal
10-16	< 10 %	10-18 %	19-23 %	> 23 %	58-72 %
17-39	< 12 %	12-20 %	21-25 %	> 25 %	53-67 %
40-55	< 13 %	13-21 %	22-26 %	>26 %	47-61 %
56-85	< 14 %	14-22 %	23-27 %	>27 %	42-56 %

For muscle mass there are no universal standard values.

## Changing the battery

A battery change is necessary when “Lo” or nothing at all is shown on the display.

When changing the battery, bear the following in mind:

- The data that has been ascertained and entered up to that point remains stored on the scale.
- Only use the type of battery that is specified in the “Technical Data”.
- Pay attention to the correct polarity (+/-) when inserting the battery. The plus pole (+) must point upwards.
- Dispose of the old battery in an environmentally responsible way (see “Disposal”).

## Cleaning

### NOTE on material damage

- Do not immerse the scale in water or other liquids.
- Do not under any circumstances use harsh, scraping, or abrasive cleaning agents or the like to clean the scale. If you do, the surface could be scratched.
- Wipe down the scale as needed using a lightly dampened, soft cloth.
- Disinfect the weighing surface and the sensors from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

## Fault/Remedy

Fault	Possible Cause
No function	<ul style="list-style-type: none"> <li>▪ Battery empty or inserted the wrong way around?</li> </ul>
The scale appear to be displaying false weighing results	<ul style="list-style-type: none"> <li>▪ Is the scale standing on a firm, even surface?</li> <li>▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scale?</li> <li>▪ Remove the battery from the scale and wait a moment. Reinsert the battery. The data stored up to this point will be retained.</li> </ul>
"Err" is shown on the display	<ul style="list-style-type: none"> <li>▪ The maximum capacity of the scale (180 kg) has been exceeded.</li> </ul>
"Lo" is shown on the display	<ul style="list-style-type: none"> <li>▪ The battery is empty (see "Changing the battery").</li> </ul>
"ErrL" is shown on the display	<ul style="list-style-type: none"> <li>▪ The fat proportion % is too low: Below 4% nothing can be displayed</li> </ul>
"ErrH" is shown on the display	<ul style="list-style-type: none"> <li>▪ The fat proportion % is too high: Above 60% nothing can be displayed.</li> </ul>
"C" is shown on the display	<ul style="list-style-type: none"> <li>▪ The weighing/measurement process was faulty: The scale are not standing stably or you are not standing still on the weighing surface during weighing.</li> </ul>

## Technical Data

---

Model:	Body Analyser Scale BA 1300-1304
Battery:	1x CR2032/3 V 
Rated current:	< 3 mA
Measuring range:	max. 180 kg
Variance:	2% +/- 100 g
Indicator intervals	
Weight and bone mass:	100 g
Body water, muscle mass, and body fat:	0.1 %
Memory cells:	12
Age input:	10-100 years
Height:	100-250 cm
Dimensions of scale:	approx. 302 x 302 x 17 mm
Distributor: (no service address)	Waagen-Schmitt GmbH Neuer Höltigbaum 15 22143 Hamburg, Germany

Our products are constantly being developed and improved.  
For this reason, design and technical modifications are possible at any time.

## Declaration of Conformity

---

Waagen-Schmitt GmbH hereby declares that the body analyser scale BA 1300-1304 complies with EU directives 2014/53/EU and 2011/65/EU.

The full text of the EU declaration of conformity is available at the following internet address: [www.ade-germany.de/doc](http://www.ade-germany.de/doc)

The declaration loses its validity if modifications have been made to the device without our approval.

Hamburg, March 2021  
Waagen-Schmitt GmbH

## Guarantee

---

The company Waagen-Schmitt GmbH guarantees for 2 years from the purchase date, the cost-free rectification of defects on the basis of material or fabrication errors through repair or exchange.

In the event of a guarantee claim, please return the scale with the proof of purchase (stipulating the reason for claim) to your dealer.

## Disposal

---

### Disposing of the packaging

Dispose of the packaging with similar materials. Put cardboard with waste paper, foil with recyclable materials.

### Disposing of the scale

Dispose of the scale in accordance with the regulations that apply in your country.

#### **Devices must not be disposed of in normal household waste.**

At the end of its life, the scale must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Hand in the old device at a collection point for electronic waste or at a recycling depot.

For more information, contact your local waste disposal firm or your local administrative body.

#### **Recyclable and disposable batteries do not belong in household waste.**

All batteries must be handed in at a collection point in your community, suburb, or at the retailer. Such batteries can thus be disposed of in an environmentally friendly way.

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