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Safety instructions for lithium batteries (CR2032/CR2430):


1. Lithium batteries are lithium manganese cells. If you do not use the scales over a longer period, you should take out the batteries of the device.
2. Never open the lithium battery, do not put into fire, or expose to shocks because an explosion risk lasts or discharge of poisonous exposure to fumes might possible.
3. Never touch leaked cells with bare hands.
4. When eyes or hands are contaminated, it is essential to rinse with lots of water. You should go to a doctor if skin or eyes are irritated.
5. Dispose cells only completely discharged or protected against polarity reversal wrapped according to your local disposal regulation.
6. Do not expose the cells to direct solar radiation and store away from heat; otherwise, there is risk of overheating.
7. Maintain a storage temperature of < 40°C.

Disposal of old electrical and electronic equipment:

This symbol on the product or its packaging indicates that this product must not be treated as normal household refuse; rather you must dispose of it at a collection point for recycling electrical and electronic equipment. You can get further information from your local council, the municipal recycling authorities or the shop where you bought the product.

Yours truly
 ADE (GmbH & Co.)

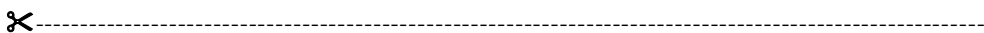Guarantee card – 3 years guarantee for the body analysing scale

<ul style="list-style-type: none"> • Sender <hr/> <hr/> <hr/>	<ul style="list-style-type: none"> • Reason for complaint <hr/> <hr/> <hr/>
--	--

• Date of purchase

Please return this guarantee card together with the receipt to your retailer.

**Model Livia / BA 1105
 Operating Instructions**

Table of Contents

1. How does the Body Analyser work?.....page 2
2. How to use.....page 2
3. Commissioning..... page 3
4. Programming / Enter of Personal Data.....page 4
5. Weighing Function (Weight Measurement only).....page 4
6. Body Analysis.....page 5
7. Error Messages..... page 6
8. Technical Data..... page 7
9. Guarantee..... page 7

BODY FAT

BODY WATER

MUSCLE MASS

BONE MASS

Thank you for purchasing this high-quality product made by ADE. You have acquired a set of body analyser scales that meet the highest technical standards and are of proven quality. Please read these operating instructions carefully before using the product for the first time and keep them in a safe place where you can find them whenever you need the information contained in the instructions.

1. How does the body analyser work?

Through these new revolutionary scales, you will be able to get information about your body composition with reference to your weight, body fat %, body water %, muscle mass % and bone mass kg. First of all, you have to enter once your personal data like gender (male or female), age and body height under your selected personal memory location. Afterwards a measurement is as simple as on a normal bathroom scale. You have only to select your personal memory location, step on the scales and within a few seconds, your measured values will be indicated. The analysis based on the measurement of the body's electrical resistance. A low, safe electrical signal passes freely through fluids contained in lean tissue, such as muscle and blood, but meets resistance passing through fat tissue. The body analyser accurately measures the resistance and with reference to your personal data use it to calculate elements of body composition.

2. How to use

If possible, the measurement should be done undressed but always barefoot.

By nature, a body is subject of fluctuations (for example, loss of water after taking exercise or after a sauna visit, after a meal, etc.). To enable comparable readings it is important to take the measurements as possible always under the same conditions.

In order to ensure that the results of analysis are accurate and consistent, please keep the measurement conditions constant always at the same daytime. In the morning, directly after you've woken up, your body will be dehydrated and also the majority of fluid will be stored in the central trunk area. This means that the resistance is relatively high because of the poor conductivity in the body. By intensify moving during the day progresses the fluid becomes more evenly distributed through the limbs and increasing the accuracy of the readings. To get a realistic result it is important to take your readings at the same time and under the same conditions every day. The easiest way to stick the rule, measure always in the morning, 15 minutes after you've got up, either before or after you've gone to the bathroom.

In spite of a very low electrical signal, this body analyser is not suitable for persons with heart pacemaker or other electronic medical implants! There is no risk for persons with none electronic implant for instance a titanium plate in the knee. However, metal conducts the electric current extremely well and therefore the scales implicates the body fat value "beautified" – means the actual body fat value is higher as indicated. However, tendencies may be determined as well.

Due to menstrual cycle and hormonal changes women experience greater changes in

8. Technical Data

Capacity x graduation:	max. 150 kg x 100 g
Display resolution:	body fat: 0.1%
Display resolution:	body water: 0.1%
Display resolution:	muscle mass: 0.1%
Display resolution:	bone mass: 100 g
Personal memory loc.:	10 positions
Age:	10 – 80 years
Body height:	100 – 220 cm
Dimensions:	335 mm x 297 mm x 25 mm
Batteries:	1 x CR2032 (included)

Subject to technical changes due to further-development of the product.

9. Guarantee

ADE guarantees to rectify any defects due to material or manufacturing defects by means of repair or replacement free of charge for a period of three years. When making a guarantee claim, please return the scales together with the guarantee card and the sales slip to your retailer and stating the reasons for the claim.

CE conformity


This device is interference-suppressed in accordance with the applicable EU directive 2004/108/EC

Note: Extreme electromagnetic influences, e.g. when operating radio communication equipment is being operated next to the scales, can affect the displayed value. When the interference is over, you can use the product for its intended purpose again; it may be necessary to switch it on again.

Statutory notice of liability for disposal of batteries:

Batteries are not domestic refuse. As a consumer, you are legally liable to return used batteries. You can return your used batteries to the public collection points in your area or anywhere where batteries of the respective type are sold.

Note:

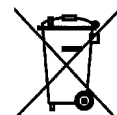You will find these markings on the batteries:

Li = The battery contains lithium

Al = The battery contains alkali

Mn = The battery contains manganese

CR (Li); A (Al, Mn); AA (Al, Mn); AAA (Al, Mn)

The average total body water percentage range for healthy men is between 50 – 65%. The equivalent value for women is approximately 45 – 60% due to the higher body fat content. About two third of the total body water is stored in the cells and is therefore called intra cellular. The other third of the total body water is extra cellular body water.

For men is the normal range of body mass index between 20 and 25, for women 19 to 24. When people in the second half of life a body mass index to 26 is considered as still normal.

Women 					
AGE	Body Fat				Body Water
	slim	normal	slightly overweight	overweight	normal
10-16	<18%	18-28%	29-35%	>35%	57-67%
17-39	<20%	20-32%	33-38%	>38%	47-57%
40-55	<23%	23-35%	36-41%	>41%	42-52%
56-85	<24%	24-36%	37-42%	>42%	37-47%

Men 					
AGE	Body Fat				Body Water
	slim	normal	slightly overweight	overweight	normal
10-16	<10%	10-18%	19-23%	>23%	58-72%
17-39	<12%	12-20%	21-25%	>25%	53-67%
40-55	<13%	13-21%	22-26%	>26%	47-61%
56-85	<14%	14-22%	23-27%	>27%	42-56%

There are no validated standards for muscle mass.

7. Error Messages

- “EEEE“ = Overload display – The 150 kg maximum capacity of the scales has been exceeded. Immediately remove the load from the scales; otherwise, the load cell could be damaged.
- “LO“ = Indication that the battery is flat – Please replace the battery with one CR2032.
- “Errl“ = Low fat indication - Under 5% body fat, not possible to indicate.
- “ErrH“ = High fat indication – Over 75% body fat, not possible to indicate.
- “.....“ = Incorrect weighing- / measuring process - For example, the scales is placed unstable or you do not remain stationary on the weighing platform.

their hydration level. Retaining fluid may cause the weight to fluctuate day-to-day. This will also affect the body fat percentage readings. If the hydration level is high, the body fat percentage is slightly less – a less hydration level means the body fat percentage is slightly higher. There is no healthy risk for pregnancy. The body analyzer is absolute safe for the mother and her unborn baby. However, due to the changes in water content in the mother’s body, the readings should not be interpreted as completely accurate.

ADE Germany takes no responsibility for any damages or losses caused by the use of this body analyzer scales or for Third Party claims. This product is exclusive intended for private use. This product is not designed for commercial or professional use in hospitals or other medical facilities.

3. Commissioning

Safety guidelines: The appliance is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless a person responsible for their safety has given them supervision or instruction concerning use of the appliance. Children should be supervised to ensure that they do not play with the appliance.

Not suitable for persons with electronic implants (heart pacemaker, etc.).

Place the scales on a plane, firm surface (tiles, hardwood flooring, etc.) to avoid a wrong measurement on carpeted floor

Attention - Hazard of slips! Do not stand on the scales with wet feet or on a slippery surface of the weighing platform.

Ensure that no object fall onto the scales, the glass surface is breakable.

Only use a damp cloth to clean the scales. Never use solvents or abrasive cleanings. Do not submerge the scales under water.

Preparation:

Open the lid of the battery compartment on the bottom of the scales and insert the supplied battery (one CR2032). Please ensure the correct polarity of the battery (+/-).

There is a little push button left beside the battery compartment on the bottom of the scales for country specific settings with the option to switch from KG = kilogram to LB = pound. Please ensure that the scales are set to the appropriate units of weight for you.

4. Programming / Enter of Personal Data

Before first time use, your personal data like gender (female or male), age, and body height must be programmed and saved under selected locations 01 – 10 (up to 10 persons). To enter your personal data use the buttons right beside the weighing platform: ON (on button), SET (save button), UP (plus button), DOWN (minus button) and follow the particular points: (numbers quick-run through continuous press):

1. To switch-on the scales press the ON-button. The display indicates “8888” and after the scales have stabilised, they display “0.0 kg”.
2. Now, press the SET-button. The display flashes USER “01”. To select your personal location (01 - 10) press the UP- or DOWN-button. After the selection of your personal location, press the SET-button to save the reading.
3. After you have selected your personal location, the symbol for gender starts to flash. To select the requested gender (♂ male or ♀ female) press the UP- or DOWN- button. After the selection of your gender, press the SET-button to save the reading.
4. After you have selected your gender, the symbol for age starts to flash (default age 25). To enter your age, press the UP- or DOWN-button. After the selection of your age, press the SET-button to save the reading.
5. After you have selected your age, the symbol for body height starts to flash (default 165 cm). To enter your body height, press the UP- or DOWN-button. After the selection of your body height, press the SET-button to save the reading.

The scales indicate once again your selected personal data. That settles the completion of input and the scales return to the indication “0.0 kg”. The scales are now ready to use and you can start your body analysis. Without further activity, the scales automatically powers-off. Your personal data are stored.

To enter the personal data of further persons repeat the complete steps listed below.

5. Weighing Function (Weight Measurement only)

1. Switch-on the scales by pressing the ON-button. The display indicates “8888” and after the scales have stabilised, they display “0.0 kg”.
2. Just, step on the scales. Remain stationary on the weighing platform and do not hold you tight.
3. Once the scales have stabilised, your weight will indicate successively three times and then for 10 seconds steady. Afterwards, the scales switches-off automatically.

6. Body Analysis

An analysis is only possible in bare feet. Your personal data must be entered first (please see page 4).

1. Switch-on the scales by pressing the SET-button. Select by using the UP- or DOWN- button the personal memory location where you stored your personal data and confirm each by pressing the SET-button. The display indicates your personal data and afterwards “0.0 kg”. The scales are now ready and you can start the body analysis.
2. Now, carefully step in bare feet on the right and left stainless steel electrodes (measuring units) of the scales and remain stationary.
3. First, your weight indicates on the display. Remain stationary on the scales. Shortly after, while the analysis is taking place, the signal “o” starts to move on the display. After completion of the measurement, your measured values will be shown successively three times as follows:
 1. your weight in kg
 2. your body fat in % and beneath an assessment by subtitle:
Underfat, Healthy, Overfat, Obese
 3. your body water in %
 4. your muscle mass in %
 5. your bone mass in kg

Afterwards, the scales switches-off automatically.

The body fat content is shown as the amount of stored fat as a percentage of total body weight. There is no generally accepted set value as this varies significantly according to gender (due to the different build) and age. Body fat is vital for basic body functions. It's protecting organs, cushioning joints, regulating body temperature, storing vitamins and serves the body as energy storage. As essential for women a minimum body fat proportion between 10% and 13% and for men between 10% and 13% is requested - with the exception of competitive athletes.

Body fat is not a clear indicator for health. Too high body weight and a too high body fat proportion are linked to the incidence of many civilization diseases such as diabetes, cardiovascular diseases, etc. Lack of movement and improper nutrition are often the prime course for these diseases. There is an obvious direct relationship between overweight and lack of movement.

A significant part of the human body consists of body water. This body water is differently distributed throughout the whole body. Lean muscle tissue contains approx. 75% body water, blood approx. 83% body water, body fat approx. 25% body water and bones have a water content of about 22%.

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