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## Declaration of conformity by the manufacturer


These scales have been manufactured in accordance with the harmonized European standards. They conform to the regulations of the EC-directives listed below:

- EMC directive 2004/108/EC regarding electromagnetic compatibility
- RoHS Directive 2011/65/EC in the currently valid versions.

This declaration loses its validity if modifications are made to the scale without our approval.

Hamburg, May 2015

**WAAGEN-SCHMITT GmbH**  
Hammer Steindamm 27 – 29  
D-22089 Hamburg
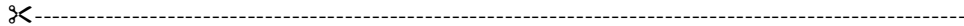
### Guarantee card – 3 years on the body analysis scales

• Sender

---

---

---

• Reason for Claim

---

---

---

• Date of purchase

Please return this section and proof of purchase to your retailer.

BODY FAT

BODY WATER

MUSCLE MASS

BONE MASS

CALORIE  
REQUIREMENT

AUTO OFF

# ADE

## Model Linda / BA 1500 Operating Instructions

### Table of Contents

1. How do body analysis scales work? .....	Page 2
2. Conditions of Use .....	Page 2
3. First Use .....	Page 3
4. Programming / Inputting personal data .....	Page 4
5. Weighing function (only weight measurement) .....	Page 4
6. Body analysis .....	Page 5
7. Error Messages .....	Page 6
8. Technical Data .....	Page 6
9. Warranty .....	Page 7

Thank you for purchasing these high-quality body analyser scales.

Please read these Operating Instructions carefully before using the product for the first time, and keep them in a safe place so that you can consult them again if you need further information.

## 1. How do body analyser scales work?

These new, revolutionary scales provide information about your body composition in terms of weight, body fat %, body water %, muscle mass %, bone mass kg and your basal metabolic rate (energy requirement) in kcal. To start, enter all your personal data, e.g. gender (male or female), height and age. This information must be stored in your personal memory location just once. After that, the process is as easy as on any regular scales. You only need to choose your individual memory button, step onto the scales and in just a few seconds your current body status will be displayed. A small electrical current is used to measure the resistance. Water, which is mainly contained in muscles, is a good conductor, fat on the other hand blocks the current. Taking the measured resistance and the stored data into consideration, the body analyser scales can very precisely calculate your body composition.

## 2. Conditions of use

The measurements should take place without clothes, if possible, and always barefoot.

Since the body is subject to natural fluctuations (e.g. dehydration through sports or sauna, intake of food and fluids), constant conditions should be ensured when measuring, so as to receive comparable values.

Care should be taken that the measurements are always carried out at the same time of day. Directly after getting up in the morning, the water has not yet spread through the body. This means that the resistance is still relatively high, due to the bad conductivity in the body. With increasing movements throughout the day, the water spreads through the body, increasing conductivity, which causes the resistance, and thus the fat level too, to be shown as lower. To determine the actual changes of body composition, it is therefore important to measure always at same time. An easy way of doing this is by always measuring a quarter of an hour after getting up in the morning, either before or after using the toilet.

**Despite the minimum current, persons with pacemakers or other implants, which emit frequencies, must not use the body analyser scales!** Other implants, e.g. a titanium plate in the knee, pose no danger. Metals, however, conduct electricity extremely well, which has the effect that body fat values are “beautified”, meaning that the actual fat level is higher than shown by the scales. Tendencies, on the other hand, can be determined just as well. In women, menstrual factors have an effect on the

## 9. Warranty


WAAGEN-SCHMITT agrees to rectify any material or manufacturing defects by offering free repair or replacement, for a period of 3 years from the purchase date. At the time of purchase, please get the retailer to fill in and stamp the guarantee card. In the case of guarantee claims, please return the scales to your retailer including guarantee card and stating the reasons for the claim.

**Note:** Extreme sources of electromagnetic radiation, such as those caused by the operation of radio equipment in the vicinity of the scales, can affect the value displayed. Once the interference has been eliminated, the product will once again function normally. Note: the unit may need to be reset.

### Legal information obligation regarding the disposal of batteries

Do not dispose of batteries in domestic waste. As a consumer, you are legally required to return used batteries. You may return your used batteries to public collection points in your community or to any retailer that sells batteries of the same type.

#### Note:

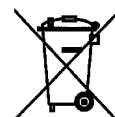The following symbols are used on batteries:

Li = Battery contains lithium

Al = Battery contains alkaline

Mn = Battery contains manganese

CR (Li); AA (Al, Mn); AAA (Al, Mn)

### Disposal of electrical and electronic devices

The symbol on the product or its packaging indicates that this product must not be disposed of with regular household waste; instead, it must be disposed of at a collection point for the recycling of electrical and electronic equipment. For further details about recycling, contact your local council, municipal recycling companies or the shop where you bought the product.

Women 					
Age	Body fat				Body water
	slim	normal	slightly overweight	overweight	normal
10-16	< 18 %	18-28 %	29-35 %	> 35 %	57-67 %
17-39	< 20 %	20-32 %	33-38 %	> 38 %	47-57 %
40-55	< 23 %	23-35 %	36-41 %	> 41 %	42-52 %
56-85	< 24 %	24-36 %	37-42 %	> 42 %	37-47 %

Men 					
Age	Body fat				Body water
	slim	normal	slightly overweight	overweight	normal
10-16	< 10 %	10-18 %	19-23 %	> 23 %	58-72 %
17-39	< 12 %	12-20 %	21-25 %	> 25 %	53-67 %
40-55	< 13 %	13-21 %	22-26 %	> 26 %	47-61 %
56-85	< 14 %	14-22 %	23-27 %	> 27 %	42-56 %

There are no general guide values for muscle mass.

## 7. Error messages

- "Err"** = Overload – The 180 kg maximum load capacity has been exceeded. Remove the object from the scales immediately to avoid damaging the scales' load cell unit.
- "Lo"** = Batteries are empty – Please replace with the appropriate battery type (3 x 1.5V AAA).
- "Err 2"** = the body-fat percentage and the body-water percentage cannot be determined and displayed, since the values are either too low (under 5%) or too high (over 60%).

## 8. Technical data

Load capacity x gradation: Max. 180 kg x 100 g  
 Display resolution: Percentage of body fat: 0.1 %  
 Display resolution: Percentage of body water: 0.1 %  
 Display resolution: Percentage of muscle mass: 0.1 %  
 Display resolution: Bone mass: 100 g  
 Pers. memory locations: 10  
 Age: 10 – 80 years  
 Height: 100 – 220 cm  
 Dimensions: 303 mm x 303 mm x 24 mm  
 Batteries: Three x 1.5V AAA (included)

Subject to technical changes.

measurement results, since hormonal fluctuations influence the body's water balance. If the body contains all lot of water, the fat values are lower, if less water is contained, the values are a little higher. When a woman is pregnant, there is no danger for the unborn child. But due to larger water deposits during pregnancy, the values are distorted making the measurement senseless.

WAAGEN-SCHMITT does not assume liability for damages or losses, caused by the body analyser scales, or for claims of third parties. This product is intended for private use only. It is not designed for professional use in hospitals or medical facilities.

## 3. First use

**Safety guidelines:** This device must not be used by individuals (including children) with limited physical, sensory or mental abilities, or who are lacking experience and understanding, unless they are supervised by another individual, who will be responsible for their safety, or who have been instructed in the use of the device. It should be ensured that children do not play with the device.

Not suitable for persons with electronic implants (cardiac pace makers, etc.).

Place the scales on a flat and solid surface (tiles, parquet, etc.). Incorrect measurements make occur on carpets.

Attention! There is a danger of slipping if the scales' surface is wet or if feet are wet.

Ensure that no objects fall onto the scales; otherwise, the glass might break.

Use only a damp cloth to clean the scales. Never use solvents or abrasive cleaning agents. Do not immerse in water.

### Preparation:

Open the battery compartment cover on the underside of the scales and insert the provided batteries. Please be mindful of the correct polarity of the batteries (+/-).

### Units of weight:

Touch the SET button as long as necessary until the kg indicator blinks on the display. By pressing the ▲- or ▼-button, you can now change the country-specific weight unit (kg = kilogram / lb = pound or st:lb = stone pound). Please make sure that the scales are set to the weight unit relevant to you.

#### 4. Programming / Inputting personal data

Prior to first use, your personal data – such as sex (female or male), age, and height – must be entered and saved in one of the memory cells P0 – P9 (for up to 10 people). The scales make use of highly sensitive sensor buttons (Sensor Touch). The touch keys are located below the display. To enter your personal data use the buttons: SET (power-on and memory button), ▲ (plus button), ▼ (minus button) and follow the individual points (quick scan by touching the button for more than 1 second):

1. Gently touch the SET button. Links im Display blinkt die Speicherplatzanzeige (z.B. P1). To choose the desired memory cell (P0 – P9), touch the ▲- or ▼- button and then save by touching the SET button.
2. After you have chosen the memory cell (e.g. memory cell P3), the symbol for sex will blink. In order to choose the desired sex (♂ male or ♀ female), touch the ▲- or ▼-button and then save by pressing the SET button.
3. Once you have chosen your sex (e.g. ♀ female), the indicator for height blinks on the display (default: 170 cm). In order to enter your height, touch the ▲- or ▼-button and then save by touching the SET button.
4. Once you have entered your height, the indicator for age blinks on the display (default age: 25) In order to enter your age, touch the ▲- or ▼- button and then save by touching the SET button.

Once you have finished entering your data, the scales return to "0.0 kg". The scales are now ready and you can undertake a body analysis. Without further activity the scales deactivate automatically after a few seconds. Your personal data has been saved.

#### 5. Weighing function (weight measurement only)

1. The scales have a step-on function. Simply step onto the scales, keep still and do not hold onto anything.
2. Once the scales have stabilized, your weight will be displayed twice in succession and then steadily for approx. 10 seconds. The scales then automatically turn off.

#### 6. Body analysis

An analysis is only possible when you are barefoot. The personal data must be entered beforehand (see section 4).

1. Activate the scales by touching the SET button. Initially, the personal data of the last setting appears on the display and a few seconds later "0.0 kg". If this is not your data, choose the memory cell in which you saved your personal data using the ▲- or ▼-button. Your personal data will be displayed and then the scales return to "0.0 kg".

2. Now stand carefully on the scales, barefoot with your feet left and right on the stainless steel strips (measuring units), and remain still.

3. First you will be shown your weight. After that, while the body analysis is being carried out, the signal "o" will run on the display. Remain standing still. After your measurements have been taken, the measured values will be shown to you twice on the display, as follows:

1. Your weight in kg
2. **BF** = your body-fat proportion in %
3. **BW** = your body water in %
4. **BON** = your bone mass in kg
5. **MUS** = your muscle mass in %
6. **KCAL** = your basal metabolic rate (daily calorie consumption)

The scales will switch off automatically after a few seconds.

You can display the values of your most recent measurements, once you have completed your initial analysis process. Simply select your assigned memory location (1-12) as described under step 1. The memory location display will stop flashing after a few seconds, and the display will show "0.0 kg". You can now press ▲ to start the automatic display of your most recent measurement values (the values are shown in sequence), while pressing ▲ or ▼ will allow you to select the values manually. Pressing the SET button will return the display to analysis mode.

The percentage of body fat shows the percentage of fat in relation to the total weight of a body. There is no general recommendation for this value, since this depends significantly on gender (due to the different physique) and age. Body fat is vital for the daily body functions. It protects organs, cushions joints, regulates body temperature, stores vitamins and acts as an energy storage for the body. With the exception of athletes, a body fat percentage for women between 10 % and 13 % and for men between 2 % and 5 % is seen as vital.

The percentage of body fat is not a clear indicator for health. A body weight and a percentage of body fat, which are too high, are related to many diseases of civilization such as diabetes, cardiovascular diseases, etc. A lack of exercise and the wrong diet are the most common causes for these diseases. There is a clear causal connection between being overweight and a lack of exercise.

A significant part of the human body consists of water. This body water has a different distribution throughout the body. Lean muscle tissue contains about 75 % water, blood approx. 83 % water, body fat approx. 25 % water and the bones have a water percentage of about 22 %. Water makes up approx. 60 % of men's body mass. The value for women is about 55 % (due to the higher percentage of body fat). Of this total amount of body water, approx. two thirds are stored in the cells and is therefore called intracellular water. The remaining third is extracellular water.

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