## Digitale Babywaage

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# Dear Customer!

You have chosen to purchase a high-quality product from the ADE brand, which combines intelligent functions with an exceptional design. This baby scale accurately records your child's weight. The many years of experience of the ADE brand ensure a high technical standard and proven quality.

We wish you many relaxing moments!

Your ADE Team

## General information

---

### About this manual

This operating manual describes how to safely use and care for your products. Keep this operating manual in a safe place in case you would like to refer to it later. If you pass this item on to someone else, pass on this manual as well.

Failure to observe the recommendations in this operating manual can lead to injuries or damage to your product.

### Explanation of symbols

This symbol, combined with the word DANGER, warns against serious injury.

This symbol, combined with the word WARNING, warns against moderate and mild injuries.

This symbol, combined with the word NOTE, warns against material damage.

This symbol signifies additional information and general advice.

# Table of Contents

---

General information .....	12
Intended Use .....	14
Safety .....	14
At a glance .....	16
Scope of delivery .....	16
Start-up .....	17
Weighing .....	18
Switching off the scales.....	18
Resetting the scales.....	18
Cleaning.....	19
Fault/Remedy.....	19
Technical Data .....	20
Declaration of Conformity.....	20
Warranty .....	20
Disposal.....	21

## Intended Use

---

The baby scale is intended for weighing babies or small animals. The scale is designed exclusively for private use and is not a medical device.

Therefore, the scale may not be used in clinics, doctor's surgeries, etc.

The maximum carrying capacity of 20 kg must not be exceeded.

## Safety

---

- This item can be used by children 8 years of age and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge, if they are supervised or have been instructed in the safe use of the item and understand the dangers that may arise.
- Children must not play with the item.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries, when swallowed, can be life-threatening. You should therefore keep the batteries and scales out of reach of small children. If a battery is swallowed, medical assistance must be sought immediately.

## DANGER for Children

- Keep children away from packaging material. Risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.

## RISK of fire/burning and/or explosion

- There is a risk of explosion if batteries are replaced incorrectly. Replace the battery using only the same or an equivalent type of battery. Pay attention to the "Technical Data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Otherwise there is an increased risk of explosion and of leakage; furthermore, gases may escape.
- Avoid contact with skin, eyes, and mucosa. If you come into contact with battery acid, flush the affected site immediately with plenty of clean water and see a doctor promptly.

## **WARNING against risk of injury**


- Do not leave the child lying, standing or sitting unattended on the scale. Hold the child firmly, if necessary. Make sure you can intervene at any time to prevent the child from slipping or falling off the scale.

## **BEWARE of material damage**

- Place the scales on a stable, even surface.
- Do not overload the scales. Pay attention to the “Technical Data”. Place yourself on the scales exclusively for weighing purposes.
- Protect the scales from sharp jolts and vibrations.
- Remove the batteries from the scales if they have been used up or if you are not going to use the scales for a long time. You will thus avoid damage that may result from the leakage of battery acid.
- When inserting and changing batteries, pay attention to the correct polarity (+/-).
- Do not subject the battery to any extreme conditions, by – for example – storing it on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Only insert batteries of the same type; do not use different types or used and new batteries with each other.
- Do not make any modifications to the item. Only allow repairs to be carried out at a specialist workshop. Repairs carried out incorrectly can considerably endanger the user.
- Do not place the scales in the vicinity of devices that may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indicators or malfunctions.
- There are feet on the underside of the scales. Since surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place an anti-slip underlay beneath the scales.

## At a glance

---

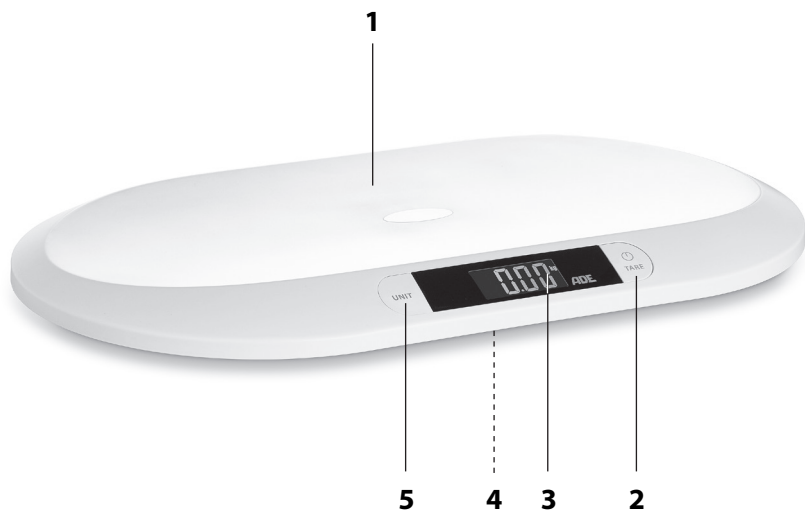- 1 Weighing surface
- 2 /TARE button:  
Switch the scales on and off, set the values to zero (tare function)
- 3 Display
- 4 Battery compartment (underside)
- 5 UNIT button: Changing the unit of measurement

## Scope of delivery

---

- Baby Scale
- Battery LR03 (AAA), 1.5 V  (2x)
- Operating Manual

# Start-up

## Inserting/changing batteries

For use, you require 2 batteries of type LR03 (AAA), 1.5 V (included). The battery compartment is located on the underside of the device.

1. Press the locking tab slightly towards the battery compartment cover and remove the battery compartment cover marked with **OPEN**.
2. Insert the batteries as shown on the base of the battery compartment. Pay attention to the correct polarity (+/-).
3. Put the battery compartment cover back in place. The locking tab must snap securely into place.

Notes on changing the battery:

- A change of batteries becomes necessary when “LO” appears on the display, or when nothing at all is displayed.
- Only use the type of battery that is specified in the “Technical Data”.
- Always change both batteries at the same time.
- Dispose of batteries in an environmentally appropriate way (see “Disposal”).

## Choosing the unit of measurement

When the scales are delivered, the unit of measurement is set to “kg” as standard. However, the measurement result can also be displayed in the English units of measurement, “lb” (pound) or “st” (stone:pound).

1. To switch on the scale, press the /TARE button.
2. Repeatedly press the UNIT button. Note the display.

## Positioning the scales

### **WARNING** against injuries

- Be sure to keep the baby safe on a firm and stable base surface of the scale. The scale and base surface, e.g. a table, must not wobble in any case.
- Place the scales on an even, firm, and dry surface. Do not place the scales on a carpet.

# Weighing

---

## WARNING against injuries

- Do not leave the child lying, standing or sitting unattended on the scale. Hold the child firmly.
- Make sure you can intervene at any time to prevent the child from slipping or falling off the scale.

1. Press the /TARE button.

The illuminated display will show "HELO" and "0.00" shortly thereafter, as well as the set unit of measurement, e.g. "kg".

2. Place a soft towel or small blanket on the weighing surface.

3. Press /TARE to reset the scale to zero (observe the display indication "TARÉ"). This ensures that the towel or blanket is not weighed.

4. Place the baby to be weighed on the towel or blanket.

The current weight value will be displayed dynamically shortly thereafter.

- Of course, you can also perform the weighing without any towel or blanket. In this case, steps 2 and 3 will be omitted.
- Weigh the baby as undressed as possible.
- Try to weigh yourself always at the same time and under consistent conditions in order to obtain comparable values.

## Switching off the scales

---

- Hold down /TARE until the indicator on the display goes out.

The scales automatically switch off if they are not used for approx. 20 seconds.

## Resetting the scales

---

Reset the scales to factory settings if unusual or obviously incorrect values are shown on the display.

1. Remove the batteries from the battery compartment and wait a moment.
2. Re-insert the batteries.

Pay attention to the correct polarity of the batteries (+/-).

After switching the scales on, you can weigh yourself as usual.

# Cleaning

## NOTE on material damage

- Do not immerse the scales in water or other liquids.
- Do not under any circumstances use harsh, scraping, or abrasive cleaning agents or the like to clean the scales. If you do, the surface could be scratched.
- Wipe down the scales as needed using a lightly dampened, soft cloth.
- Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

## Fault/Remedy

Fault	Possible Cause
No function	<ul style="list-style-type: none"> <li>▪ Battery empty or inserted the wrong way around?</li> </ul>
The scales appear to be displaying false weighing results	<ul style="list-style-type: none"> <li>▪ Are the scales standing on a firm, even surface?</li> <li>▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scales?</li> <li>▪ Remove the battery from the scales and wait a moment. Reinsert the battery.</li> </ul>
"O-Ld" is shown on the display	<ul style="list-style-type: none"> <li>▪ The maximum capacity of the scales (20 kg) has been exceeded.</li> </ul>
"LO" is shown on the display	<ul style="list-style-type: none"> <li>▪ The batteries are empty (see "Inserting/changing batteries").</li> </ul>
"UnSt" is shown on the display	<ul style="list-style-type: none"> <li>▪ The scale was unstable during the measurement. Check the base of the scale and repeat the weighing.</li> </ul>
"- - -" is shown on the display	<ul style="list-style-type: none"> <li>▪ The scale is not ready or is positioned incorrectly. Reposition the scale.</li> </ul>

## Technical Data

---

Model:	Digital baby scale BE 1817
Battery:	2x LR03 (AAA), 1.5 V 
Rated current:	20 mA
Measuring range:	max. 20 kg
Accuracy:	1 % +/- 10g
Indicator division:	10 g
Dimensions of scales:	approx. 557 x 320 x 35 mm
Distributor: (no service address)	Waagen-Schmitt GmbH Hammer Steindamm 27-29 22089 Hamburg, Germany

Our products are constantly being developed and improved.  
For this reason, design and technical modifications are possible at any time.

## Declaration of Conformity

---

Waagen-Schmitt GmbH hereby declares that digital baby scale BE 1817 complies with EU directives 2014/30/EU and 2011/65/EU.

The declaration loses its validity if modifications have been made to the device without our approval.

Hamburg, February 2018  
Waagen-Schmitt GmbH  
Hammer Steindamm 27-29  
22089 Hamburg, Germany

## Warranty

---

The company Waagen-Schmitt GmbH guarantees for 2 years from the purchase date the cost-free rectification of defects on the basis of material or fabrication errors through repair or exchange.

In the event of a warranty claim, please return the product with the proof of purchase (stipulating the reason for claim) to your dealer.

# Disposal

---

## Disposing of the packaging

Dispose of the packaging with similar materials. Put cardboard with waste paper, films with recyclable materials.

## Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

### **Devices must not be disposed of in normal household waste**

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Hand in the old device at a collection point for electronic waste or at a recycling depot.

For more information, contact your local waste disposal firm or your local administrative body.

### **Recyclable and disposable batteries do not belong in household waste**

All batteries must be handed in at a collection point in your community, suburb, or at the retailer. Such batteries can thus be disposed of in an environmentally friendly way.

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