Personenwaage BE 1613/1614 **FITvigo**

Bathroom Scales | Báscula personal
Pèse-personne | Bilancia pesapersona
Waga osobowa

FITvigo

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Dear Customer!

You have chosen to purchase a product from the FITvigo program from ADE: Activity trackers, bathroom and body analyser scales, as well as blood pressure monitors are equipped with a Bluetooth function in order to transfer all data to the FITvigo app. FITvigo is free and, alongside clear progress monitoring, offers additional motivation and stimulation for more wellbeing. At the same time, the many years of experience of the brand ADE guarantees the high quality of the product and the app.

We hope you enjoy your FITvigo product.

Your ADE Team

General Information

About this manual

This operating manual describes how to safely use and care for your products. Keep this operating manual in a safe place in case you would like to refer to it later. If you pass this item on to someone else, pass on this manual as well.

Failure to observe the recommendations in this operating manual can lead to injuries or damage to your product.

Explanation of symbols

This symbol, combined with the word DANGER, warns against serious injury.

This symbol, combined with the word WARNING, warns against moderate and mild injuries.

This symbol, combined with the word NOTE, warns against material damage.

This symbol signifies additional information and general advice.

Table of Contents

General Information	14
Intended Use.....	16
Safety	16
At a glance (scope of delivery)	18
Scope of delivery	18
Start-up	19
Weighing	20
Use FITvigo app.....	21
Cleaning.....	21
Fault/Remedy.....	22
Technical Data	23
Declaration of Conformity	23
Guarantee	24
Disposal	24

Intended Use

The bathroom scales are designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

The maximum carrying capacity of 180 kg must not be exceeded.

Safety

- This item can be used by children 8 years of age and above and by persons with reduced physical, sensory, or mental abilities or lacking experience and knowledge, if they are supervised or have been instructed in the safe use of the item and understand the dangers which may arise.
- Children must not play with the item.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries, when swallowed, can be life-threatening. You should therefore keep the batteries and scales out of reach of small children. If a battery is swallowed, medical assistance must be sought immediately.

DANGER for children

- Keep children away from packaging material. Risk of choking if swallowed. Also beware that children do not pull the packaging material over their heads.

RISK of fire/burning and/or explosion

- Risk of explosion if batteries are improperly exchanged. Replace the batteries only with one of the same or a similar type. Pay attention to the "Technical Data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited.
- Avoid contact with skin, eyes, and mucosa. If contact is made with the battery fluid, rinse the affected areas immediately with plenty of clean water and see a doctor straight away.


WARNING against risk of injury

- Do not step onto the outer edge of any one side of the scales. Risk of tipping!
- Do not step onto the scales with wet feet or when its surface is damp. Risk of slipping!
- The weighing surface is made of glass. If you fall down, there is a risk of the glass breaking and thus a risk of injury from glass splinters.

BEWARE of material damage


- Place the scales on a stable, even surface.
- Do not overload the scales. Pay attention to the "Technical Data". Place yourself on the scales exclusively for weighing purposes.
- Protect the scales from sharp jolts and vibrations.
- Remove the batteries from the scales if they have been used up or if you are not going to use the scales for a long time. You will thus avoid damage that may result from the leakage of battery acid.
- Pay attention to the correct polarity (+/-) when exchanging the battery.
- Do not expose the batteries to any extreme conditions, e.g. by storing them on heaters or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the item. Only allow repairs to be carried out at a specialist workshop. Repairs carried out incorrectly can considerably endanger the user.
- There are feet on the underside of the scales. Since floor surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place an anti-slip underlay beneath the scales.

At a glance (scope of delivery)

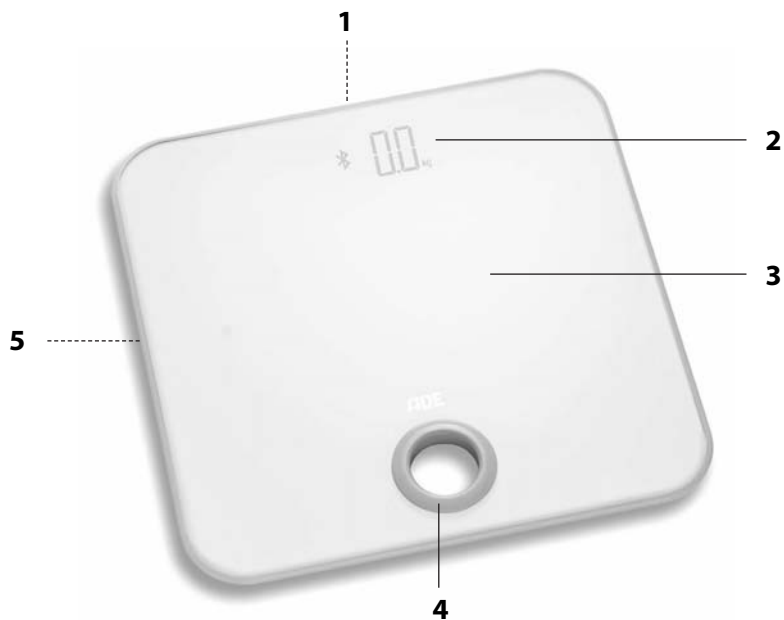- 1 Unit and Bluetooth button: Choose unit of measurement, connect scales with the app
- 2 Display
- 3 Weighing surface
- 4 Grip
- 5 Battery compartment

Scope of delivery

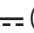- Scales
- Battery type R03 (AAA), 1.5 V  (4x)
- Operating Manual

Start-up

Inserting/changing batteries

For use you require 4 batteries of type R03 (AAA), 1.5 V (included). The battery compartment is located on the underside of the device.

1. Slide the battery compartment cover in the direction of the grip until it stops.
2. Remove the battery compartment cover.
3. Insert the batteries as depicted on the bottom of the battery compartment.
Pay attention to the correct polarity (+/-).
"8888" appears briefly on the display, then "0.0 kg".
4. Put the battery compartment cover back in place.

Notes on changing the battery:

- A change of batteries becomes necessary when "Lo" appears on the display, or when nothing at all is displayed.
- Only use the type of battery that is specified in the "Technical Data".
- Always exchange all of the batteries at the same time.
- Dispose of batteries in an environmentally appropriate way (see "Disposal").

Choosing the unit of measurement

When the scales are delivered, the unit of measurement is set to "kg" as standard. However, the measurement result can also be displayed in the English units of measurement, "lb" (pound) or "st:lb" (stone:pound).

- Press the Unit and Bluetooth button on the rear side of the scales repeatedly. Pay attention to the display indicator.

Install app

1. Make sure your smartphone is connected with the internet.
2. Download the "FITvigo" app in the app store (iOS) or Google Play store (Android).

- Depending on the settings of your smartphone, additional costs may be incurred for the internet connection or data transfer.
- We recommend downloading over a wireless network.

Connect the device with the app (pairing)

1. Activate Bluetooth on your smartphone.
2. Start the app and follow the instructions. Pay attention to the information and explanations displayed.
After a short period, you'll have the scales connected with the FITvigo app.

Positioning the scales

- Place the scales on an even, firm, and dry surface. Do not place the scales on carpet.

Weighing

WARNING against injuries

- Do not step onto the scales with wet feet or when its surface is damp.
Risk of slipping!

- Weighing is also possible without the FITvigo app connected.
- If you have connected the FITvigo app, the measured weight will be transmitted to the FITvigo app after the weighing has concluded.

1. With a little force, step once briefly onto the weighing surface.
"0.0" and the set unit of measurement, e.g. "kg", will be shown on the display.
2. Position yourself on the scales. Remain standing still and do not hold on to anything.
Shortly thereafter, the measured weight is displayed and the data is transferred to the FITvigo app.
A few seconds after you step off the weighing surface, the scales will switch off automatically.

When the scales are switched off you can also step onto the scales directly.

Use FITvigo app

The data ascertained by the scales are transmitted to the FITvigo app, as long as you have connected it to the scales (see "Connect Device with the App").

In this regard, pay attention to the following symbols:

Symbol	Meaning
	Existing Bluetooth connection
	Data is being transmitted <ul style="list-style-type: none"> ▪ Data transmission successful: the symbol disappears. ▪ Data transmission unsuccessful: the symbol continues to be displayed. The pending measurement data is temporarily stored in the scales and, after conclusion of the next measurement, transferred to your smartphone.

Abide by the following points in order to prevent possible transmission problems:

- The distance between the scales and your smartphone should be between 1 m and 10 m.
- There should be no obstructions between the scales and the smartphone.
- Other electronic devices (in particular those that are Bluetooth enabled) should be at least 1 m from the scales.
- Ensure that the Bluetooth function on your smartphone is switched on.

In the app your calculated BMI (body mass index) is also displayed alongside your body weight.

Further information on the FITvigo app can be found online at www.fitvigo.de

Do you need help or do you have additional questions?

Then please send us an email at: support@fitvigo.de

Cleaning

NOTE on material damage

- Do not immerse the scales in water or other liquids.
- Do not under any circumstances use harsh, scraping, or abrasive cleaning agents or the like to clean the scales. If you do, the surface could be scratched.
- Wipe down the scales as needed using a lightly dampened, soft cloth.
- Disinfect the weighing surface and the sensors from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

Fault/Remedy

Fault	Possible Cause
No function	<ul style="list-style-type: none"> ▪ Are the batteries empty or have they been inserted incorrectly?
The scales appear to be displaying false weighing results.	<ul style="list-style-type: none"> ▪ Are the scales standing on a firm, even surface? ▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scales? ▪ Remove the batteries from the scales and wait a moment. Re-insert the batteries. The data stored up to this point will be retained.
“----” is shown on the display.	<ul style="list-style-type: none"> ▪ The maximum capacity of the scales (180 kg) has been exceeded.
“Lo” is shown on the display	<ul style="list-style-type: none"> ▪ The batteries are empty (see “Inserting/ changing batteries”).
While attempting to connect the device with the app (pairing), “E I” is shown on the display.	<ul style="list-style-type: none"> ▪ Has Bluetooth been activated on your smartphone? ▪ Is the FITvigo app active?
Data transmission has failed.	<ul style="list-style-type: none"> ▪ Was the distance between the scales and the smartphone too great? Max. distance: 10 m.

Technical Data

Model:	Bathroom Scales BE 1613/1614
Battery:	4x R03 (AAA)/1.5 V 
Rated current:	≤ 60 mA
Measuring range:	max. 180 kg
Variance:	0-50 kg = +/- 0.3 kg 51-100 kg = +/- 0.4 kg 101-150 kg = +/- 0.5 kg 151-180 kg = +/- 0.7 kg
Dimensions of scales:	approx. 320 x 320 x 25 mm
Distributor (no service address):	Waagen-Schmitt GmbH Hammer Steindamm 27-29 22089 Hamburg, Germany

Our products are constantly being developed and improved.
For this reason, design and technical modifications are possible at any time.

Declaration of Conformity

These scales have been manufactured in accordance with the harmonised European standards. It conforms to the provisions of the EC directives cited below:

- 2014/30/EU EMC Directive
 - 2011/65/EU RoHS Directive
- in the currently valid versions.

This declaration loses its validity if modifications are made to the product without our approval.

Hamburg, August 2016
Waagen-Schmitt GmbH
Hammer Steindamm 27-29
22089 Hamburg, Germany

Guarantee

The company Waagen-Schmitt GmbH guarantees for 3 years from the purchase date the cost-free rectification of defects on the basis of material or fabrication errors through repair or exchange.

In the event of a guarantee claim, please return the product with the proof of purchase (stipulating the reason for claim) to your dealer.

Disposal

Disposing of the packaging

Dispose of the packaging with similar materials. Put cardboard with waste paper, films with recyclable materials.

Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

Devices must not be disposed of in normal household waste.

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Hand in the old device at a collection point for electronic waste or at a recycling depot.

For more information, contact your local waste disposal firm or your local administrative body.

Recyclable and disposable batteries do not belong in household waste.

All batteries must be handed in at a collection point in your community, suburb, or at the retailer. Such batteries can thus be disposed of in an environmentally friendly way.

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